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

Distributor (US):
pro3dure medical LLC
9825 Valley View Road · Eden Prairie, MN 55344
Phone: 952-426-1928

info@pro3dure.com · www.pro3dure.com · Made in Germany

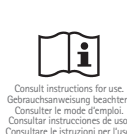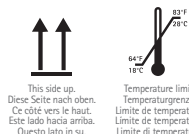Rx only
QTY: 1EA

1. Intended use / Indication

Polyurethane material for the production of earmolds.

2. Contraindication

- SysTherm[®] ECO is contraindicated ...
1. ... if a patient is known to be allergic to any of the ingredients.
 2. ... for direct contact with uncured or partially cured material.
 3. ... for any application that is not part of the indication (see above).

3. Patient Target Group

Persons for whom an earmold is to be created.

4. Intended Users

Earmold laboratory, ENT specialist, hearing aid acoustician

5. Requirements

Scale required. Store material at room temperature.

6. Material

SysTherm[®] ECO is a 2-component polyurethane-based material consisting of a base and a catalyst.

7. Note

The batch number and best-before date are indicated on each package. In case of complaints, please always state the batch number of the product. Do not use the product after the expiration date. Deviations from the described manufacturing processes or storage conditions may lead to changes in the mechanical and optical properties of the material. Due to its chemical base, SysTherm[®] ECO is sensitive to moisture and should be stored in a dry place at room temperature. In accordance with the EU Medical Devices Regulation, users/patients are obliged to report serious incidents involving a medical device to the manufacturer and the competent authority of the country in which they occurred.

8. Application

Production process based on 3D-printed molds or silicone molds using SysTherm[®] SPT-1 or SPT-2 separator; see the instructions for use for SysTherm[®] SPT-1 or SPT-2. The compatible cast materials can be found in the product compatibility table (see right).

- 1) Open the bottles (components A and B).
- 2) Pour the two components A and B in a ratio of 5:1 into a plastic cup placed on a scale. Ensure that weighing tolerances do not exceed 1%. If a tolerance is unavoidable, it is better to have a small excess of component A than of component B.
- 3) Mix the components evenly at room temperature using a suitable spatula in a plastic beaker for 4 minutes. Ensure that the material on the walls and bottom of the cup is also mixed homogeneously and that no dead spaces are created. Avoid introducing air bubbles as much as possible. If air bubbles remain in the mixture, this is not a problem as they will disappear during pressure polymerization.
- 4a) Carefully pour the mixed SysTherm[®] ECO into a silicone mold or ...
- 4b) ... into a prepared 3D-printed cast mold (see SysTherm[®] cast 2.0). Avoid air pockets and also fill the inlet as a reservoir for possible polymerization shrinkage.
- 5) Place the filled cast molds in a dry pressure pot. Make sure that SysTherm[®] ECO does not flow out. Cure SysTherm[®] ECO for 1.5 hours at 80 °C / 3 bar. Take the heating time into account.

- 6) Carefully remove the SysTherm[®] earmold from the silicone or cast mold. Wear safety goggles and use suitable tools. Clean and process the SysTherm[®] ECO earmolds. Avoid aggressive solvents such as acetone or isopropanol.

9. Post-processing

Remove the raw earmold from the silicone or 3D-printed cast mold. It can then be finished with a carbide burr. The earmold can then be sanded with sandpaper or processed using a vibratory finishing process.

10. Cleaning

Clean the finished earmold with soapy water and then dry it with a paper towel.

11. Cleaning and care

After wearing your earmold, rinse it with warm soapy water and clean all surfaces with a soft paper towel. Store the SysTherm[®] ECO earmold in a dry box for hearing aids.

12. Cleaning agents and conditions to avoid

- toothpaste or commercially available denture cleaning agents
- alcohol-based cleaning agents or mouthwashes
- abrasive cleaners or general household cleaners
- hot water (> 60 °C) or ultrasonic cleaning devices

13. Component A

Hazard Statements:

May cause an allergic skin reaction. Harmful if inhaled. May cause respiratory irritation. Toxic to aquatic life with long lasting effects.

Precautionary Statements:

Avoid breathing dust/fume/gas/mist/vapours/spray. Use only outdoors or in a well-ventilated area. Contaminated work clothing should not be allowed out of the workplace. Avoid release to the environment. Wear protective gloves/protective clothing/eye protection/face protection. IF ON SKIN: Wash with plenty of soap and water. IF INHALED: Remove person to fresh air and keep comfortable for breathing. Get medical advice/attention if you feel unwell. If skin irritation or rash occurs: Get medical advice/attention. Take off contaminated clothing and wash it before reuse. Collect spillage. Store in a well-ventilated place. Keep container tightly closed. Store locked up. Dispose of contents/container according to official regulations.

14. Component B


Hazard Statements:

Causes serious eye irritation. Suspected of damaging fertility or the unborn child.

Precautionary Statements:

Obtain special instructions before use. Do not handle until all safety precautions have been read and understood. Avoid contact during pregnancy/while nursing. Wash hands thoroughly after handling. Wear protective gloves/protective clothing/eye protection/face protection. IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. IF exposed or concerned: Get medical advice/attention. If eye irritation persists: Get medical advice/attention. Store locked up. Dispose of contents/container according to official regulations.

UMDNS 16-697

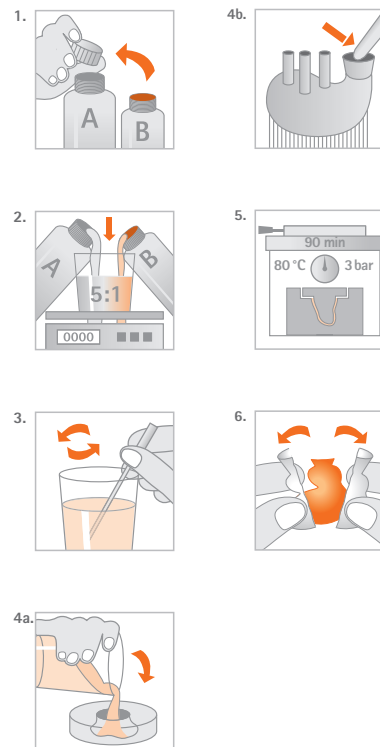
Produkt-Kompatibilität

	SysTherm [®] ECO
audioprint [®] GR-12.1	SysTherm [®] SPT-2
audioprint [®] GR-12.2	SysTherm [®] SPT-2
audioprint [®] GR-12.3	SysTherm [®] SPT-2
SysTherm [®] cast 2.0	SysTherm [®] SPT-1 or SPT-2

Physical properties*:

SysTherm[®] ECO

- Setting time:
90 min. at
80 °C / 176 °F
in pressure pot
at 3 bar
- Final hardness:
Shore A
at 23 °C: 80-95
Shore A
at 37 °C: 65-75
- Application:
at 23 ± 2 °C,
50 ± 5 %
relative humidity
- Storage:

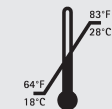
Ordering information:

SysTherm[®] ECO

- A-component**
1000 g bottle
REF: A4020250
- B-component**
clear-transparent
200 g bottle
REF: A4020252
- rose-transparent**
200 g Flasche
REF: A4020254
- A-component**
500 g bottle
REF: A4020251
- B-component**
clear-transparent
100 g bottle
REF: A4020253
- rose-transparent**
100 g bottle
REF: A4020255

* These data come from measurements of a representative sample that was determined as part of our quality assurance.

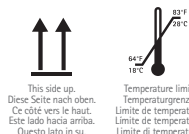Do not use if package is damaged.
Nicht verwenden, wenn das Paket beschädigt ist.
Ne pas utiliser si l'emballage est endommagé.
No lo use si el paquete está dañado.
Non usare se il pacchetto è danneggiato.

Consult instructions for use.
Gebrauchsanweisung beachten.
Consulter le mode d'emploi.
Consultar instrucciones de uso.
Consultare le istruzioni per l'uso.

Keep away from sunlight.
Von Sonnenlicht fernhalten.
Éloigner du soleil.
Mantener alejado de la luz solar.
Tenere lontano dalla luce del sole.

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pro3dure medical GmbH